

ECG 300GT


Digital Three Channel ECG

Connected with PC In many ways

Provide USB interfaces to have the real-time transmission, with the of computer and going for the data remote diagnosis and analysis management and print through the PC






Output illustrations of Diagnostic report

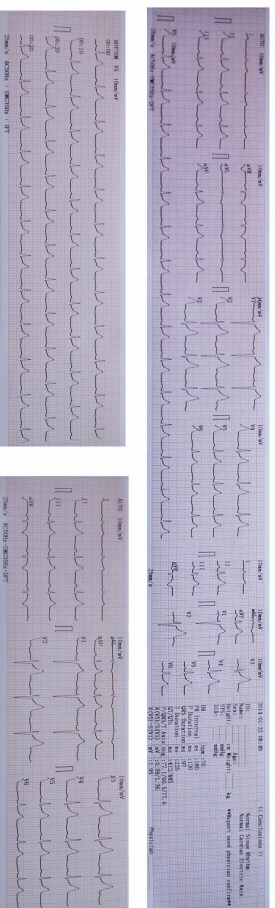

Touch Color LCD Screen



Accessories

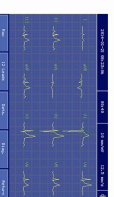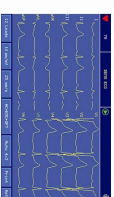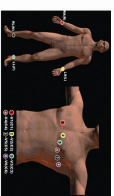

Features

- 4.3" color TFT touch screen displays ECG waveform and work status clearly, achieves the function of observing firstly, then printing, which saves the recording paper.
- Sync collection for 12-lead ECG, digital signal processor, to get high-quality ECG waveforms via AC filter, baseline filter and EMG filter of ECG signal.
- Display 3/6/12-lead ECG, print mode, sensitivity, paper speed, filter states in a screen, convenient for contrastively diagnosing.
- Print modes and formats are optional, such as manual, auto 4×3, auto 3×4+1, auto 3×4, auto 3-2+1, auto 3-2, auto 2×6+1, auto 2×6, auto 1×12-1, auto 1×12, rhythm 4, rhythm 3, rhythm 2, etc; time for print waveform can be adjusted, and timing print function meets different application requirements.
- With the functions of auto-analysis and auto-diagnosis for routine ECG parameters, provide measurement results and auto-diagnosis conclusion for HR, P-R interval, P Duration, QRS Duration, Q-T interval, Q-Tc, P Axis, QRS Axis, T Axis, RV5, SV1), RV5+SV1), etc, which reduces the doctor's burden.
- A/C/D/C, in optional DC state, built-in rechargeable lithium battery can continuously work for 10-hour, continuously print more than 3-hour and up to 260-ECG, which meets the requirements of examination, visiting a patient at home.
- Built-in memorizer can store up to 1000 cases, which is convenient for doctor to review case and statistic.
- Optional operation interface with different languages: Chinese, English, Spanish, Turkish, Polish, Italian, French, German, Portuguese, Russian, which meets different requirements.

ECG 300GT



CONTEC MEDICAL SYSTEMS CO.,LTD.
<http://www.contecm.com>